


ALTERNATING PRESSURE ON-DEMAND LOW AIR LOSS

OEM Series



MATTRESS SYSTEM

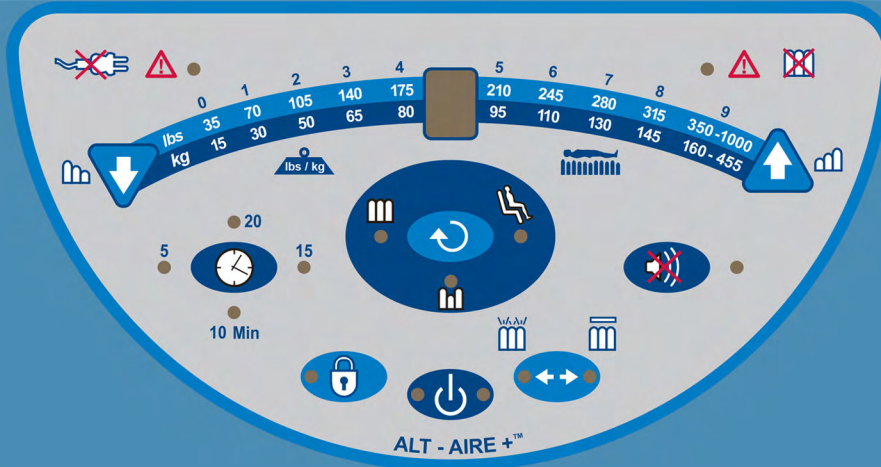


FOAM AIRE SYSTEM



BOLSTER SYSTEM



MADE IN THE USA

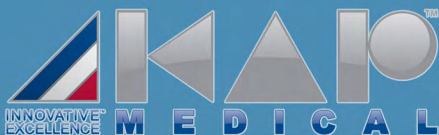

K-2oem ALT-AIRE +™

K-2oem ALT-AIRE +™

Alternating Pressure with On-Demand Low Air Loss (LAL) System

• PATENTED • Smart Support Surface™ with Real Time Pressure Monitoring (RTPM) sensor technology

Effective means of providing ALTERNATING (AP) low interface pressure therapy with On-Demand Low Air Loss relief. Simple to operate and affordable. Smart Support Surface™ (SSS) provides real time and accurate patient interface pressures below capillary occlusion with enhanced Alternating Pressure - Low Air Loss therapy and patient comfort.

Control Unit (K-2oem) OEM

- Plug and play auto-sensing multi-use control unit, operates in any country. Single high airflow pump, 42~52 LPM (liters per minute). State of the art micro-controller with Real Time Pressure Monitoring (RTPM) sensor technology. Software flexibility allows customer changes (custom pressures, times, and functions).
- User-friendly control panel with dome keys and LED display for simple user interface. Displays current mode, A/P times, and patient comfort pressure level from Soft 0 to Firm 9 (7 ± 5 to 32 ± 5 mmHg).
- Micro-controller shuts off the pump once the desired patient comfort pressure is obtained. The pump is off most of the time. Max Flow mode with 30 minute timer for rapid inflation of the mattress, inflates in 5 ~ 25 minutes depending on the size of the mattress.
- **ON-DEMAND LAL:** Low Air Loss relief can be turned on or off. During on demand LAL mode the system provides continuous LAL relief.
- Provides effective static or alternating pressure therapies along with continuous low air loss relief. In upright mode the pump inflates the mattress to 80% higher pressure to eliminate patient bottoming. **Optional Wireless Auto Upright (Fowler) Sensor (K-150).** Solenoids in the pump hold air for days in the event of power failure or during patient transportation.
- Lock out key completely locks the control panel to eliminate tampering of the patient settings, including power switch.
- Audiovisual alarms in the event of power outage or low pressures (hose disconnect). Alarm silence key to mute alarms. Internal memory always retains current settings during power outage, will resume desired operation upon power restoration.
- Control unit is reliable, quiet, compact and lightweight. For quick inflate please use our K-39 RID™ Rapid Inflator / Deflator unit.


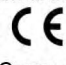
Smart Support Surface™ (K-2oemM) OEM, Patented Low Air Loss Top Sheet


- Therapeutic support surface consists of superior quality, urethane coated nylon air cushions (16 cell overlay or 16/20 cell mattress) to provide effective A/P therapy with low interface pressures. Pad overlays or cell overlays are 3" or 5" or 7" high and mattress replacements are 8" or 10" high.
- **OPTIONAL SIDE BOLSTERS** constant air pressure side bolsters / foam bolsters help contain the patient on the mattress. Air bolsters can be deflated for easy patient access into and out of the mattress and also for easy patient procedures.
- **OPTIONAL PATENTED LOW AIR LOSS** 3 layered zippered top sheet. Top layer is low friction and low shear force producing, and highly vapor permeable urethane coated breathable nylon or 4-way stretch material, the middle layer is a low air loss spacer fabric layer and the bottom layer is a liquid proof layer. 2" convoluted safety foam base or 2" optional air pad base (8" & 10" mattress) to provide additional support to the patient in the event of power failure.
- The air cells are connected to the manifold using quick connectors (hose manifold with straight connectors on Foam Aire products), and fastened to the base. The top sheet is zippered to the base. Tri-Lumen tubing assembly with single, rugged quick disconnect connector for rapid attach and CPR deflate function. Holds air for days during power failure.
- **FOAM AIRE**, K-2oemFAM mattress has fire barrier Kevlar lining with multi-ply visco foam top and hard perimeter bolster which houses the foam filled air cells.
- Complete support surface is made of superior quality materials and is modular in construction. All of the components such as manifold, tubing assembly, air cushions, top sheet, and foam base are interchangeable and can be easily cleaned or detached for laundry.
- Lightweight support surface can be rolled up easily and has two handles or two shoulder straps, which can be used to carry by hand or over the shoulder.

SAC™ Support Air Cushion (K-140)

- SAC™ (optional) foot support air cushion can be used to support patient's feet and to provide zero interface pressures in the heel area.

ALT-AIRE +™ MODEL #: K-2oemCOS / K-2oemMS8 / K-2oemMS / K-2oemRSB / K-2oemFAMS TECHNICAL DATA

AUTO- SENSING.	K-2oem CONTROL UNIT	SUPPORT SURFACE K-2oemCO / K-2oemM8 / K-2oemM / K-2oemFAM / K-2oemRSM
Input Voltage	90~264 V, 50/60 Hz, 1A/0.5A, Universal, AUTO-SENSING	System Model # K-2oemMS for 80"x36"x10" Mattress replacement system
Power Cord	14' (427cm) long hospital grade detachable.	3"/5"/7" (7.6/12.7/17.8cm) high mattress overlay.
Fuse Protection	1 Amp Slow blow dual fuses.	8"/10" (20.5/25cm) high mattresses. Side bolsters are available
A/P Cycle Time	5, 10, 15, and 20 minutes.	36"/39"/42"/48"/54"/60" (91/99/107/122/137/152cm) wide mattresses.
Dimensions (LWH) / Weight	13.5" x 11" x 5.5" (34cm x 28cm x 14cm) / 13 lbs. (5.8 Kg.).	76"/80"/84" (193/203/213cm) long mattresses.
Weight Capacity	Standard & Bariatric 360 ~ 1000 lbs (163~455 Kg).	Standard Size: 80" x 36" x 5" (Overlay Mattress), 80"x36"x10" (Mattress Replacement) / 23 lbs. (10.5 Kg.). Custom sizes available.
Safety Agency Approvals	 ETL LISTED CONFORMS TO IEC STD 60601-1 IEC STD 60601-1-2 EN STD 60601-1 ANSI/AAMI STD ES60601-1 CERTIFIED TO CAN/CSA STD C22.2 NO. 60601-1 Intertek 2003311	With respect to Electrical Shock, Fire and Mechanical Hazards. CE Mark: 3 rd Edition  MADE IN USA An FDA Registered Company, Product Is FDA Listed. Medicare & Medicaid coded, Health Canada Licensed. ISO 13485 certified company.
Warranty	One year limited, parts and labor.	One year limited, parts and labor.

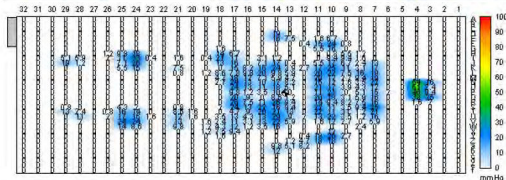

Call us or visit our website for interface pressure map.

KAP Medical products are patented and/or covered by pending patents. Trademarks referenced are property of KAP Medical. © 2013, 1395 Pico Street, Corona, California 92881 • USA
 J Within the USA 1-866-KAP-MED 1 (1-866-527-6331) • Outside the USA 1-951-340-4360 • www.kapmedical.com • sales@kapmedical.com • P/N: 400245 Rev: A ECO 13-0352